
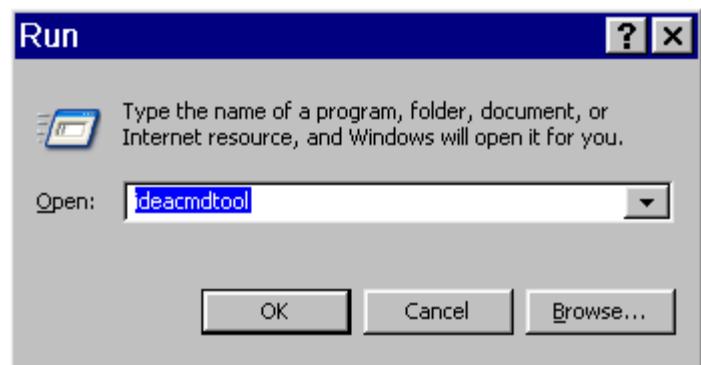


Signal Logging VD11A:

Remark: all inputs need to be confirmed by the *Enter* (CR) Key. The resulting logfiles can be found in the directory c:\MedCom\log. The extension of the ECG files is .ecg. Start *ideacmdtool*:

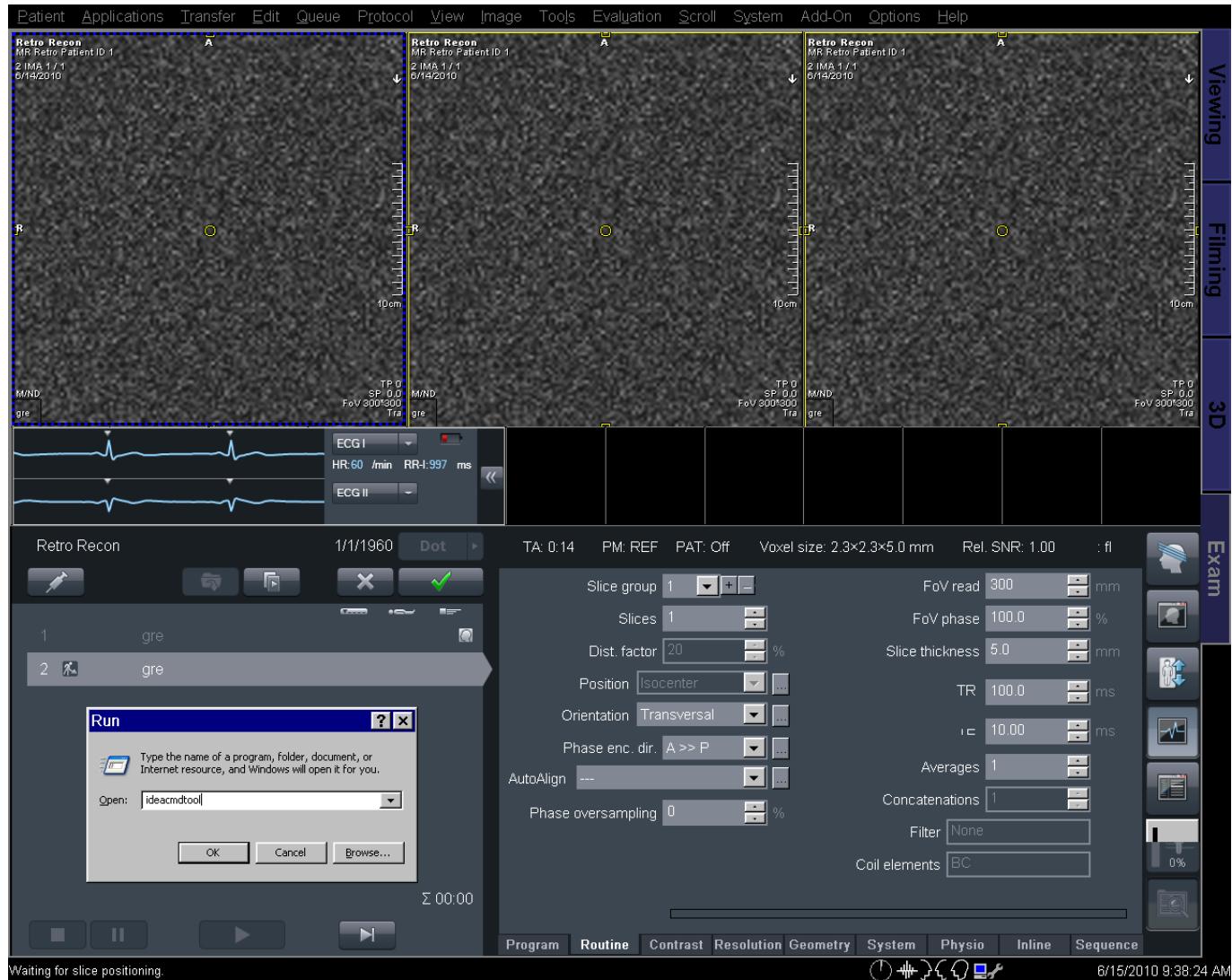

IDEACmdTool Menu

1. PMU control...
2. PAS unload

Signal Logging VD11A:

- 3. Measure
- 4. Online Export Defaults
- 5. Switches...
- 6. Set Default

Please enter number or unique text of a cmd. ('q' for quit)

> 1

PMU control Menu

- 1. PMU Signal simulation
- 2. PMU Signal logging

Please enter number or unique text of a cmd. ('q' for quit)

PMU> 2

Menu

- 1. dump
- 2. setHWTraceLevel
- 3. isCallAllowed
- 4. configureIF
- 5. startHostLogECG
- 6. stopHostLogECG
- 7. startHostLogRespiration
- 8. stopHostLogRespiration
- 9. startHostLogPulse
- 10. stopHostLogPulse
- 11. startHostLogExternal
- 12. stopHostLogExternal
- 13. startHostLogAll
- 14. stopHostLogAll

Please enter number or unique text of a cmd. ('q' for quit)

Cmd > 5

enter filename: patient1

command startHostLogECG returned: kPerRetOK

Menu

- 1. dump
- 2. setHWTraceLevel
- 3. isCallAllowed
- 4. configureIF
- 5. startHostLogECG
- 6. stopHostLogECG
- 7. startHostLogRespiration
- 8. stopHostLogRespiration
- 9. startHostLogPulse
- 10. stopHostLogPulse
- 11. startHostLogExternal

Signal Logging VD11A:

- 12. stopHostLogExternal
 - 13. startHostLogAll
 - 14. stopHostLogAll
-

Please enter number or unique text of a cmd. ('q' for quit)

Cmd > 6

command stopHostLogECG returned: kPerRetOK

Menu

- 1. dump
 - 2. setHWTraceLevel
 - 3. isCallAllowed
 - 4. configrelF
 - 5. startHostLogECG
 - 6. stopHostLogECG
 - 7. startHostLogRespiration
 - 8. stopHostLogRespiration
 - 9. startHostLogPulse
 - 10. stopHostLogPulse
 - 11. startHostLogExternal
 - 12. stopHostLogExternal
 - 13. startHostLogAll
 - 14. stopHostLogAll
-

Please enter number or unique text of a cmd. ('q' for quit)

Cmd > q

PMU control Menu

- 1. PMU Signal simulation
 - 2. PMU Signal logging
-

Please enter number or unique text of a cmd. ('q' for quit)

PMU> q

IDEACmdTool Menu

- 1. PMU control...
 - 2. PAS unload
 - 3. Measure
 - 4. Online Export Defaults
 - 5. Switches...
 - 6. Set Default
-

Please enter number or unique text of a cmd. ('q' for quit)

> q

Signal Logging VD11A:

ideacmdtool exits...